

COMPONENTS

Kit Component	Amount
96-well plate precoated with anti-rat Activin A antibody	1 Plate
Protein Standard: Lyophilized recombinant rat Activin A	2 tubes, 10 ng/tube
Sample Diluent Buffer	30 ml
Biotinylated Antibody (Anti-rat Activin A)	130 µl (100x)
Antibody Diluent Buffer	12ml
Avidin-Biotin-Peroxidase Complex (ABC) Solution	130 µl (100x)
ABC Diluent Buffer	12 ml
Tetramethylbenzidine (TMB) Color Developing Agent	10 ml
TMB Stop Solution	10 ml


Washing Buffer (not provided): TBS or PBS

0.01M TBS: Add 1.2g Tris, 8.5g NaCl; 450µl of purified acetic acid or 700µl of concentrated hydrochloric acid to 900ml H₂O and adjust pH to 7.2-7.6. Adjust the total volume to 1L.

0.01M PBS: Add 8.5g sodium chloride, 1.4g Na₂HPO₄ and 0.2g NaH₂PO₄ to 900ml distilled H₂O and adjust pH to 7.2-7.6. Adjust the total volume to 1L.

Storage

Store at 4°C. Cell Applications, Inc. recommends using the kit within 6 months of order.

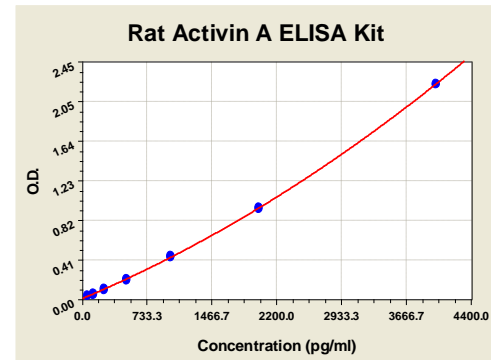

X	pg/ml	0.0	62.5	125	250	500	1000	2000	4000
Y	O.D.450	0.020	0.046	0.065	0.113	0.213	0.451	0.952	2.233

Figure 1: Activin A Standard Curve. Using the rat Activin A ELISA Kit, O.D. data was graphed against Activin A protein concentration. The TMB reaction was incubated at 37°C for 37 min.

BACKGROUND

Activin A is a homodimer of 14kDa beta-A. Activin A, a cytokine member of the transforming growth factor-beta superfamily, is expressed locally by the mesenchymal component of the hemopoietic microenvironment. Its expression is regulated on the mRNA level by different cytokines, and the biological activity of the protein is tightly controlled by several inhibitory molecules.¹ Inhibins and activins are members of the transforming growth factor beta superfamily and are known to modulate the growth and differentiation of several cell types.² Inhibins and activins inhibit and activate, respectively, the secretion of follitropin by the pituitary gland. Inhibins/activins are involved in regulating a number of diverse functions such as hypothalamic and pituitary hormone secretion, gonadal hormone secretion, germ cell development and maturation, erythroid differentiation, insulin secretion, nerve cell survival, embryonic axial development or bone growth, depending on their subunit composition. Inhibins appear to oppose the functions of activins. The standard product used in this kit is recombinant Activin A, which is composed of two single chains of 116 amino acids with the molecular mass of 26kDa.

References

1. Shav-Tal, Y. et al: Stem Cells. 20(6):493-500, 2002.
2. Di, Loreto, C. et al: Eur J Endocrinol.141(2):190-4, 1999.

ELISA OVERVIEW


Cell Applications ELISA Kits are based on standard sandwich enzyme-linked immunosorbent assay technology. Freshly prepared standards, samples, and solutions are recommended for best results.

1. Prepare test samples.
2. Prepare a protein standard of the target protein.
3. Add test samples and standards to the pre-coated 96-well plate. Do not wash.
4. Add biotinylated detection antibodies. Wash.
5. Add Avidin-Biotin-Peroxidase Complex (ABC) Solution. Wash.
6. Add Tetramethylbenzidine (TMB) Color Developing Agent, containing HRP substrate.
7. Add TMB Stop Solution
8. Subject the plate to analysis.

NOTES:

- Before using the kit, quick spin tubes to bring down all solution to the bottom of tube.
- Duplicate assay wells are recommended for both standard and sample testing.
- Do not let the 96-well plate dry, this will lead to inactivation of plate components.
- When diluting samples and reagents, ensure that they are mixed completely and evenly.
- Pre-warm diluted ABC and TMB solutions at 37°C for 30 min before use to avoid variable temperature effects.
- For washes, use TBS or PBS. Do not touch well walls.
- A protein standard is included in the kit. A protein standard detection curve should be generated with each experiment, no more than 2 hours prior to the experiment.
- The user will determine sample dilution fold by estimation of target protein amount in samples.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.

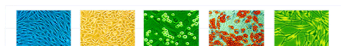

PROTOCOL

I. Plate Washing

Discard the solution in the plate without touching the side walls. Blot the plate onto paper towels or other absorbent material. Soak each well with at least 0.3 ml PBS or TBS buffer for 1–2 minutes. Repeat this process two additional times for a total of three washes.

II. Preparation of Test Samples

Test Sample Processing

- **Cell culture supernate, tissue lysate or body fluids:** Remove particulates by centrifugation.
- **Serum:** Allow the serum to clot in a serum separator tube (about 2 hours) at room temperature. Centrifuge at approximately 1000 X g for 10 min.

Sample Dilution Guideline

Estimate the concentration of the target protein in the sample and select a proper dilution factor such that the diluted target protein concentration falls within the standard curve range. Depending on the sample, several trial dilutions may be necessary. Dilute the sample using the provided diluent buffer, mixing well. Suggested working dilutions of samples are as follows:

Target Protein Concentration Range	Sample Working Dilution	Sample Vol.	Diluent Buffer Vol.
40-400 ng/ml	1:100	1 µl	99 µl
4-40 ng/ml	1:10	10 µl	90 µl
62.5-4000 pg/ml	1:2	50 µl	50 µl
≤62.5 pg/ml	n/a	100µl	n/a

If samples will be assayed within 24 hours, store at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.

III. Preparation of Reagents

Reconstitution of the Standard

The standard solutions should be prepared no more than 2 hours prior to the experiment. Two tubes of the standard are included in each kit. Use one tube for each experiment.

1. 10,000pg/ml of rat Activin A standard solution: Add 1 ml sample diluent buffer into one tube, keep the tube at room temperature for 10 min and mix thoroughly. The 10 ng/ml standard solution may be stored at 4°C for up to 12 hours, or at -20°C for up to 48 hours. Avoid repeated freeze-thaw cycles.
2. 4000pg/ml of rat Activin A standard solution: Add 0.4 ml of the above 10ng/ml Activin A standard solution into 0.6 ml sample diluent buffer and mix thoroughly.
3. 2000pg/ml→62.5pg/ml of rat Activin A standard solutions: Label 6 Eppendorf tubes with 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml, respectively. Aliquot 0.3 ml of the sample diluent buffer into each tube. Add 0.3 ml of the above 4000pg/ml Activin A standard solution into 1st tube and mix. Transfer 0.3 ml from 1st tube to 2nd tube and mix. Transfer 0.3 ml from 2nd tube to 3rd tube and mix, and so on.

Preparation of Biotinylated Antibody Working Solution

The solution should be prepared no more than 2 hours prior to the experiment.

1. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
2. Biotinylated antibody should be diluted in 1:100 with the antibody diluent buffer and mixed thoroughly.

Preparation of the Avidin-Biotin-Peroxidase Complex (ABC) Working Solution

The solution should be prepared no more than 1 hour prior to the experiment.

1. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
2. Avidin- Biotin-Peroxidase Complex (ABC) should be diluted in 1:100 with the ABC dilution buffer and mixed thoroughly.

IV. ELISA

The ABC working solution and TMB color developing agent must be kept warm at 37°C for 30 min before use. When diluting samples and reagents, they must be mixed completely and evenly. A standard detection curve should be prepared for each experiment. The user will decide sample dilution fold by crude estimation of target protein amount in samples.

1. Aliquot 0.1ml per well of the 4000pg/ml, 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml rat Activin A standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Blank well). Add 0.1ml of each properly diluted sample of rat sera, plasma, body fluids, tissue lysates or cell culture supernatants to each empty well. See “Sample Dilution Guideline” for details. We recommend that each rat Activin A standard solution and each sample is measured in duplicate.
2. Seal the plate with the cover and incubate at 37°C for 90 min.
3. Remove the cover, discard plate content, and blot the plate onto paper towels or other absorbent material. Do NOT let the wells completely dry at any time.
4. Add 0.1ml of biotinylated anti-rat Activin A antibody working solution into each well and incubate the plate at 37°C for 60 min.
5. Wash the plate three times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material.
6. Add 0.1ml of prepared ABC working solution into each well and incubate the plate at 37°C for 30 min.
7. Wash plate 5 times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1-2 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material.
8. Add 90 µl of prepared TMB color developing agent into each well and incubate plate at 37°C for 25-30 min (shades of blue can be seen in the wells with the four most concentrated rat Activin A standard solutions; the other wells show no obvious color).
9. Add 0.1ml of prepared TMB stop solution into each well. The color changes into yellow immediately.
10. Read the O.D. absorbance at 450nm in a microplate reader within 30 min after adding the stop solution.

V. Calculating Protein Concentration

- For all wells, determine O.D.450(Relative):

$$\text{O.D.450(Relative)} = \text{O.D.450(Reading)} - \text{O.D.450(Blank)}$$
- Plot the standard curve:
 Plot O.D.450(Relative) of each standard solution (Y) vs. the respective concentration of the standard solution (X). See **Figure 1** for a typical standard curve.
- The target protein concentration in samples can be interpolated from the standard curve. Multiply the interpolated concentration by the dilution factor to obtain the target protein concentration in the sample.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.

