with any other cytokine.

#### **COMPONENTS**

APPLICATIONS, INC.

Kit Component	Amount
96-well plate precoated with anti-human ANG antibody	1 Plate
Protein Standard: Lyophilized recombinant human ANG	2 tubes, 10 ng/tube
Sample Diluent Buffer	30 ml
Biotinylated Antibody (Anti-human ANG)	130 µl (100x)
Antibody Diluent Buffer	12ml
Avidin-Biotin-Peroxidase Complex (ABC) Solution	130 µl (100x)
ABC Diluent Buffer	12 ml
Tetramethyllbenzidine (TMB) Color Developing Agent	10 ml
TMB Stop Solution	10 ml

#### Washing Buffer (not provided): TBS or PBS

0.01M TBS: Add 1.2g Tris, 8.5g NaCl; 450µl of purified acetic acid or 700µl of concentrated hydrochloric acid to 900ml H<sub>2</sub>O and adjust pH to 7.2-7.6. Adjust the total volume to 1L.

0.01M PBS: Add 8.5g sodium chloride, 1.4g Na<sub>2</sub>HPO<sub>4</sub> and 0.2g NaH<sub>2</sub>PO<sub>4</sub> to 900ml distilled H<sub>2</sub>O and adjust pH to 7.2-7.6. Adjust the total volume to 1L.

#### Storage

Store at 4°C. Cell Applications, Inc. recommends using the kit within 6 months of order.

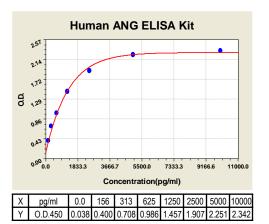


Figure 1: ANG Standard Curve. Using the Human ANG ELISA Kit, O.D. data was graphed against ANG protein concentration. The TMB reaction was incubated at 37°C for 30 min.

# **BACKGROUND**

Angiogenin (ANG)/ degranulation inhibitory protein (DIP) is a basic heparin binding protein of 14.4 kDa and belongs to a protein family, designated RNAse superfamily. It was also called ribonuclease A family 5a (Rnase5a). It has been shown that angiogenin plays an important role in tumor angiogenesis. Its expression increases in a variety of cancer cells, which results in an elevation of plasma angiogenin in cancer patients. Angiogenin is also expressed in a wide spectrum of normal cells. Angiogenin is a liverderived component of normal serum the concentration of which can increase in various disease states and its expression is regulated in vivo in a manner that is characteristic of acute phase proteins. IL6, a major inducer of acute phase proteins, stimulates the synthesis and secretion of angiogenin protein in human HepG2 cells.1

Angiogenin is now recognized as a pleiotropic molecule capable of inducing several intra- and extracellular activities. It induces most of the individual events in the process of angiogenesis. It was shown that angiogenin binds to the endothelial cell surface, interacts with a 170 kDa putative receptor or a 42 kDa binding protein on the cell surface. Angiogenin undergoes nuclear translocation in endothelial cells, which is essential for its biologic activity. It was shown that angiogenin is able to bind to the ribosomal DNA (rDNA) and to stimulate rRNA transcription. In addition, Angiogenin displays limited ribonuclease activity. The activity of angiogenin can be blocked by ribonuclease inhibitors, including a specific inhibitor of mammalians, designated RAI (ribonuclease/angiogenin inhibitor). Blocking of the ribonuclease activity concomitantly also blocks the angiogenic activity of angiogenin. The activity of angiogenin as an angiogenesis factor is blocked by Angiostatin.<sup>2</sup> In addition, human angiogenin also directly stimulates prostate cancer cell proliferation.3 Furthermore, it was shown that the serum angiogenin level increases in patients with nasopharyngeal carcinoma, and is correlate with primary tumor progression. Moreover, angiogenin has also been reported to suppress significantly the proliferation of human lymphocytes stimulated by a mixed-lymphocyte culture or by phytohemagglutinin or concanavalin A. It also inhibits both spontaneous and peptide-stimulated degranulation of polymorphonuclear leukocytes.4

## References

- Shapiro, R. et al: Biochem. 25:3527–32, 1986 Kong, H.L. & Crystal, R.G.: J. NCI 90:274-86, 1998 Yoshioka, N. et al: Proc. Natl. Acad. Sci. USA 103:14519-24, 2006
- Tschesche, H. et al: J. Biol. Chem. 269:30274-80, 1994

#### **ELISA OVERVIEW**

Cell Applications ELISA Kits are based on standard sandwich enzyme-linked immunosorbent assay technology. Freshly prepared standards, samples, and solutions are recommended for best results.

- Prepare test samples.
- Prepare a protein standard of the target protein.
- 3. Add test samples and standards to the pre-coated 96-well plate. Do not
- Add biotinylated detection antibodies. Wash.
- 5. Add Avidin-Biotin-Peroxidase Complex (ABC) Solution. Wash.
- Add Tetramethyllbenzidine (TMB) Color Developing Agent, containing 6. HRP substrate.
- Add TMB Stop Solution 7
- Subject the plate to analysis.

#### NOTES:

- Before using the kit, quick spin tubes to bring down all solution to the bottom of tube.
- Duplicate assay wells are recommended for both standard and sample testing.
- Do not let the 96-well plate dry, this will lead to inactivation of plate components.
- When diluting samples and reagents, ensure that they are mixed completely and evenly.
- Pre-warm diluted ABC and TMB solutions at 37°C for 30 min before use to avoid variable temperature effects.
- For washes, use TBS or PBS. Do not touch well walls.
- A protein standard is included in the kit. A protein standard detection curve should be generated with each experiment, no more than 2 hours prior to the experiment.
- The user will determine sample dilution fold by estimation of target protein amount in samples.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.











Target Protein Species Range Specificity Human 78pg/ml→5000pg/ml No detectable cross-reactivity with any other cytokine.

### PROTOCOL

## I. Plate Washing

APPLICATIONS, INC.

Discard the solution in the plate without touching the side walls. Blot the plate onto paper towels or other absorbent material. Soak each well with at least 0.3 ml PBS or TBS buffer for 1~2 minutes. Repeat this process two additional times for a total of three washes.

## **II. Preparation of Test Samples**

## **Test Sample Processing**

- Cell culture supernate, tissue lysate or body fluids: Remove particulates by centrifugation.
- Serum: Allow the serum to clot in a serum separator tube (about 2 hours) at room temperature. Centrifuge at approximately 1000 X g for 10 min.
- Plasma: Collect plasma using heparin, EDTA as an anticoagulant. Centrifuge for 30 min at 1000 x g within 30 min of collection. Analyze immediately or aliquot and store frozen at -20°C. Citrate is not recommended as the anticoagulant.

## Sample Dilution Guideline

Estimate the concentration of the target protein in the sample and select a proper dilution factor such that the diluted target protein concentration falls within the standard curve range. Depending on the sample, several trial dilutions may be necessary. Dilute the sample using the provided diluent buffer, mixing well. Suggested working dilutions of samples are as follows:

Target Protein	Sample	Sample	Diluent
Concentration Range	Working Dilution	Vol.	Buffer Vol.
100-1,000 ng/ml	1:100	1 µl	99 µl
10-100 ng/ml	1:10	10 µl	90 µl
156-100 pg/ml	1:2	50 µl	50 µl
≤156 pg/ml	n/a	100µl	n/a

If samples will be assayed within 24 hours, store at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.

## III. Preparation of Reagents

## Reconstitution of the Standard

The standard solutions should be prepared no more than 2 hours prior to the experiment. Two tubes of the standard are included in each kit. Use one tube for each experiment.

- 10,000pg/ml of human ANG standard solution: Add 0.5 ml sample diluent buffer into one tube, keep the tube at room temperature for 10 min and mix thoroughly.
- 2. 5000pg/ml→156pg/ml of human ANG standard solutions: Label 6 Eppendorf tubes with 5000pg/ml, 2500pg/ml, 1250pg/ml, 625pg/ml, 313pg/ml, 156pg/ml, respectively. Aliquot 0.3 ml of the sample diluent buffer into each tube. Add 0.3 ml of the above 10, 000pg/ml ANG standard solution into 1st tube and mix. Transfer 0.3 ml from 1st tube to 2nd tube and mix. Transfer 0.3 ml from 2nd tube to 3rd tube and mix, and so on.

#### **Preparation of Biotinylated Antibody Working Solution**

The solution should be prepared no more than 2 hours prior to the experiment.

- 1. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
- Biotinylated antibody should be diluted in 1:100 with the antibody diluent buffer and mixed thoroughly.

# Preparation of the Avidin-Biotin-Peroxidase Complex (ABC) Working Solution

The solution should be prepared no more than 1 hour prior to the experiment.

- 1. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
- Avidin- Biotin-Peroxidase Complex (ABC) should be diluted in 1:100 with the ABC dilution buffer and mixed thoroughly.

#### IV. ELISA

The ABC working solution and TMB color developing agent must be kept warm at 37°C for 30 min before use. When diluting samples and reagents, they must be mixed completely and evenly. A standard detection curve should be prepared for each experiment. The user will decide sample dilution fold by crude estimation of target protein amount in samples.

- . Aliquot 0.1ml per well of the 10,000pg/ml, 5000pg/ml, 2500pg/ml, 1250pg/ml, 625pg/ml, 313pg/ml, 156pg/ml human ANG standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Blank well). Add 0.1ml of each properly diluted sample of human sera, plasma, body fluids, tissue lysates or cell culture supernatants to each empty well. See "Sample Dilution Guideline" for details. We recommend that each human ANG standard solution and each sample is measured in duplicate.
- 2. Seal the plate with the cover and incubate at 37°C for 90 min.
- Remove the cover, discard plate content, and blot the plate onto paper towels or other absorbent material. Do NOT let the wells completely dry at any time.
- Add 0.1ml of biotinylated anti-human ANG antibody working solution into each well and incubate the plate at 37°C for 60 min.
- 5. Wash plate 3 times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material. (Plate Washing Method: Discard the solution in the plate without touching the side walls. Blot the plate onto paper towels or other absorbent material. Soak each well with at least 0.3 ml PBS or TBS buffer for 1~2 minutes. Repeat this process two additional times for a total of THREE washes. Note: For automated washing, aspirate all wells and wash THREE times with PBS or TBS buffer, overfilling wells with PBS or TBS buffer. Blot the plate onto paper towels or other absorbent material.)
- Add 0.1ml of prepared ABC working solution into each well and incubate the plate at 37°C for 30 min.
- Wash plate 5 times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1-2 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material. (See Step 5 for plate washing method).
- 8. Add 90 µl of prepared TMB color developing agent into each well and incubate plate at 25°C in dark for 25-30min (Note: For reference only, the optimal incubation time should be determined by end user. And the shades of blue can be seen in the wells with the four most concentrated human ANG standard solutions; the other wells show no obvious color).
- Add 0.1ml of prepared TMB stop solution into each well. The color changes into yellow immediately.
- Read the O.D. absorbance at 450nm in a microplate reader within 30 min after adding the stop solution.

## V. Calculating Protein Concentration

- For all wells, determine O.D.450(Relative):
  - O.D.450(Relative) = O.D.450(Reading) O.D.450(Blank)
- Plot the standard curve:
  - Plot O.D.450(Relative) of each standard solution (Y) vs. the respective concentration of the standard solution (X). See **Figure 1** for a typical standard curve.
- The target protein concentration in samples can be interpolated from the standard curve. Multiply the interpolated concentration by the dilution factor to obtain the target protein concentration in the sample.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.









